

# SIEMENS

## MODULARIS Uro

**CS**

### Update Instructions

**SP001/03/S**

Title: Reinforcement of tabletop

Reason for update: ☒ SafetyUrgency: ☒ Immediate ☐ Within 3 monthsUpdate materials required? ☒ Yes ☐ NoMaterials free of charge? ☒ Yes ☐ NoReturn of parts? ☐ Yes ☒ No

Estimated completion time: 0.5 h Number of CSE's: 1

Customer application training? ☐ Yes ☒ No**Systems/Products affected/System identifying IVK**

Name	Material No.	Serial No.
MODULARIS Uro	55 31 012	See Page 2 Systems/Products Affected
Uro Praktis	56 66 016	01001 - 01054

Remark: n.a.

**Components affected/to be modified/IVK**

Name	Material No.	Serial No.	Component status affected
n.a.	n.a.	n.a.	n.a.

Remark: n.a.

Chg. Ref. No.: 118 838  
Name: Erlwein  
Dept.: CS PS 24

© Siemens AG 2003

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Print No.: SPL1-130.896.03.01.02

Doc. Gen. Date: 01.03

Replaces: n.a.

Version 2.0; August 1, 2002

Page 1 of 6

**Document Revision Level**

This document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

**Disclaimer**

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such service.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

**Systems/Products Affected**

MODULARIS Uro Serial No. 01001 to 01003, 01005 to 01074, 01076 to s 01109, 01111 to 01119, 01121 to 01126.

Uro Praktis Serial No. 01001 - 01054

**Reason for the Update**

A heavier load placed on the examination area can cause the affected part of the tabletop to break, which can lead to personal injury, under certain circumstances to serious personal injury.

**Prerequisites**

n.a.

**Special Tools / Documents**

n.a.

## Ordering Information

The following update kit has to be ordered from CSML (SAP Distribution Channel, factory 2050) or BU Logistics (factory BU) (only with an order via BU Logistics):

### **Update Kit 56 67 337 (Material Number)**

List the system Serial Number on the order!

## Contents of the Update Kit

Update kit **56 67 337** contains the following parts:

Pos.	Quan.	Material No.	Name
1	1	n.a.	Support strut with screws and washers
2	1	n.a.	Stick-on label
3	1	n.a.	Technical Document

Tab. 1

## Return of Parts

n.a.

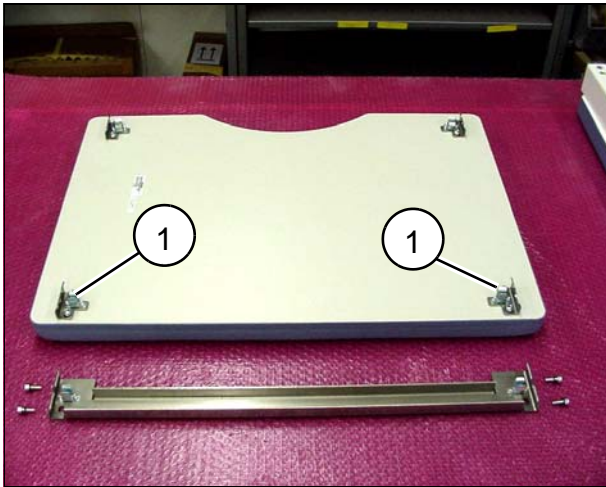


Fig. 1

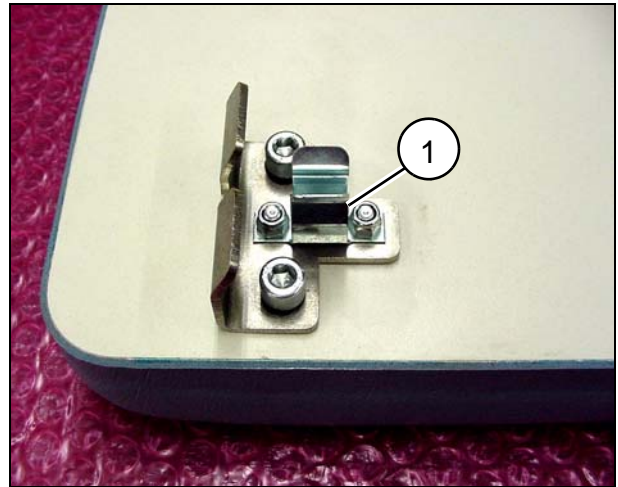


Fig. 2



Fig. 3

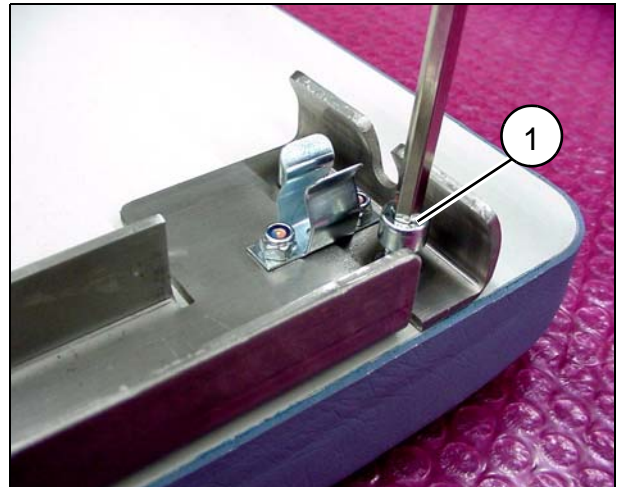


Fig. 4

## Work Steps

- Move the tabletop upwards.
- Remove both mounting clips (1/Fig. 1 and 1/Fig. 2) from the back side of the tabletop.
- Position the rail according to (Fig. 3) and loosely screw it on with the supplied screws (1/Fig. 4).
- Position the top on the table so that the clips lock in place and the strut is in alignment.
- Tighten the screws from underneath (1/Fig. 4).

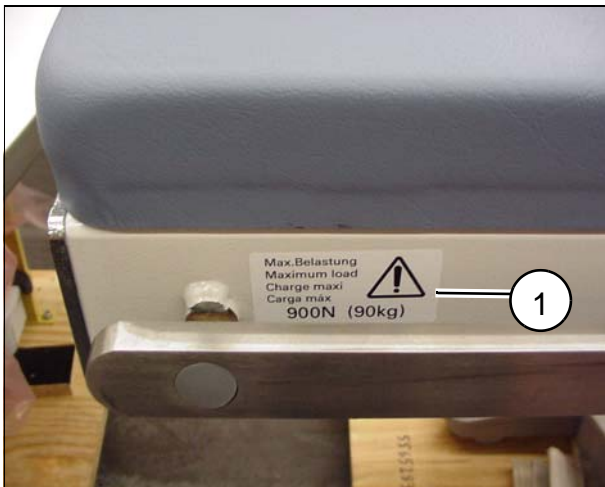


Fig. 5

- Attach the stick-on label (1/Fig. 5) on both sides of the table.

## Final Check

n.a.

## Customer Information

Inform the customer about the new load limit. The part will be replaced at a later time.

## Final Work Steps

- Update the system documentation.  
Update the revision level, the operating instructions and the technical documentation.  
Fill out, and if needed, make a copy of the attached "Completion Protocol/ Update Completion Form" and file it in the corresponding System Binder/User Handbook.
- Updates that have already been completed prior to publication of this SI must also be reported.
- The update is reported as follows:
  - The modification reply cards (Type 606) previously distributed with the publication of updates no longer apply.
  - The modification reply report has to be prepared by authorized personnel using an application on the Intranet.

## Changes to Previous Version

n.a.

## Completion Protocol

The update with the number **SP001/03/S** has been completed.

Material Number:.....

Serial Number: .....

Customer: .....

Site: .....

Customer No.: .....

Name ( CSE ): .....

Telephone: .....

Country: .....

Location: .....

Date: .....

Signature: .....

Remark: .....

.....

<b>NOTE</b>
-------------

---

**After completing the update, make a copy of this page, fill it out and file it in the corresponding System Binder/User Handbook.**

---